

TOF


DAQ test & Resin test

박관형

DAQ

- FADC is working
- Gain calibration (with 호빈, 병찬)





Optical Resin test

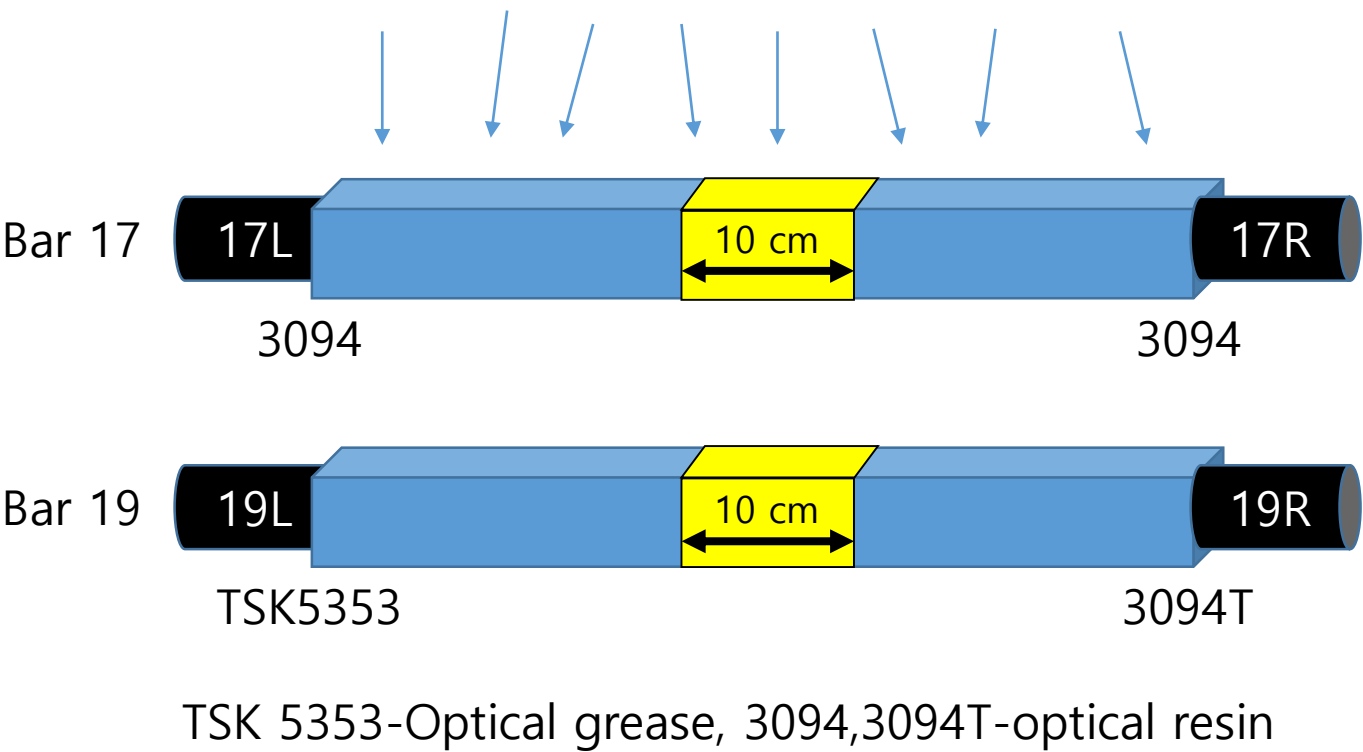


Voltage:1550 V

Date	17L	17R	19R	19L
Optical Mat.	3094	3094	3094T	TSK
18/01/03			761	699
18/01/10	1174	809		
18/01/16	1180	816		
18/01/23	1152	805		
18/01/30	1188	820		
19/02/22	570	813	581	580

Bar 17 L (optical resin)


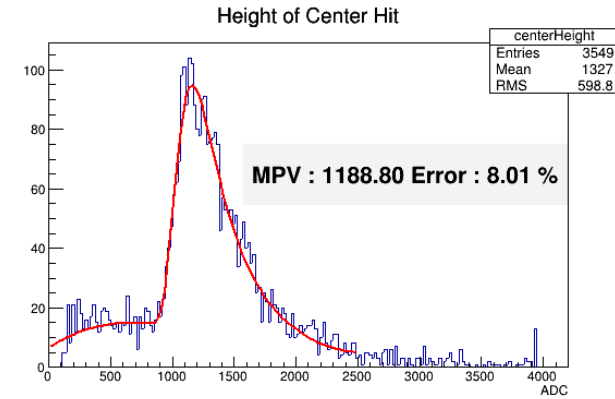
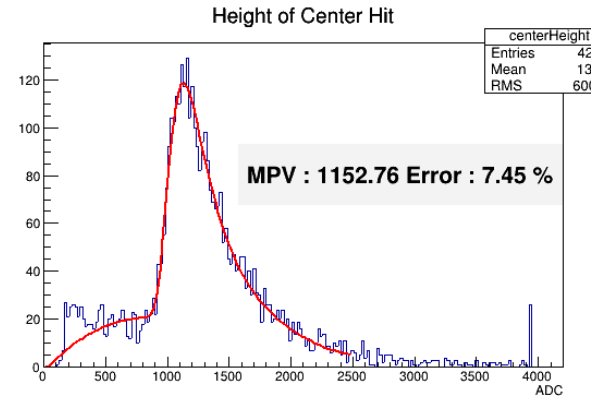
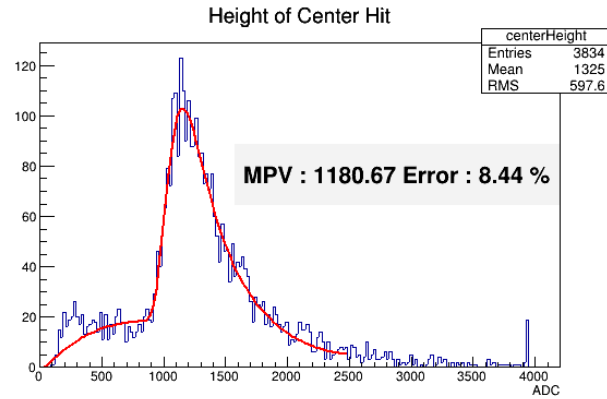
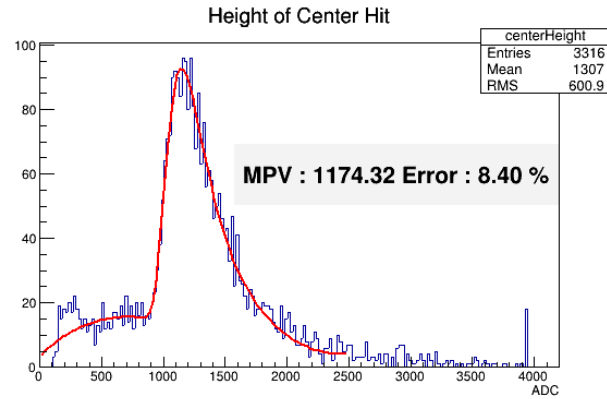


h00

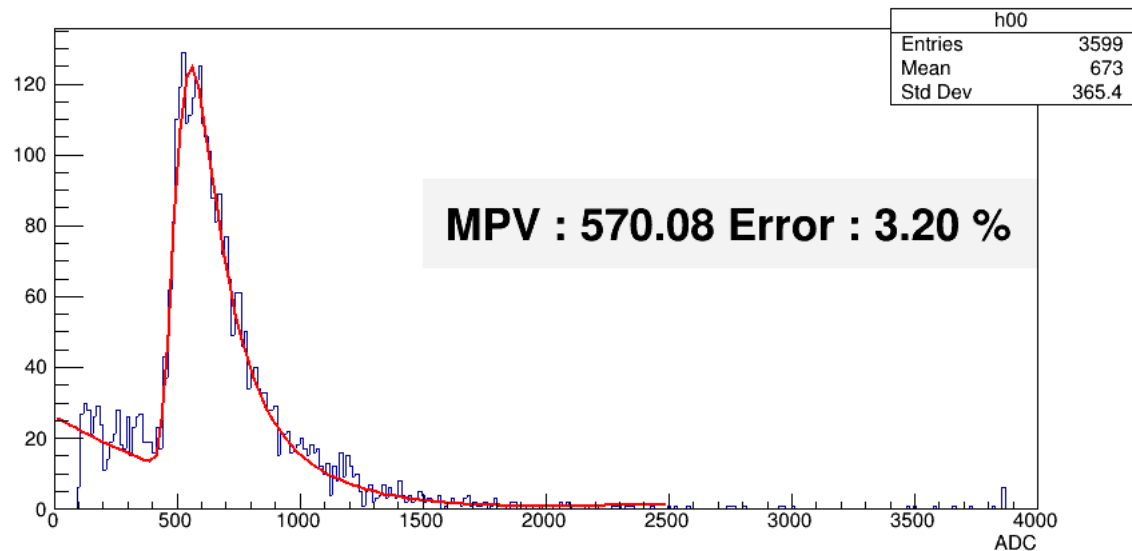



Date	MPV (Count)
18/01/10	1174
18/01/16	1180
18/01/23	1152
18/01/30	1188
19/02/22	570

Bar 17 R (optical resin)

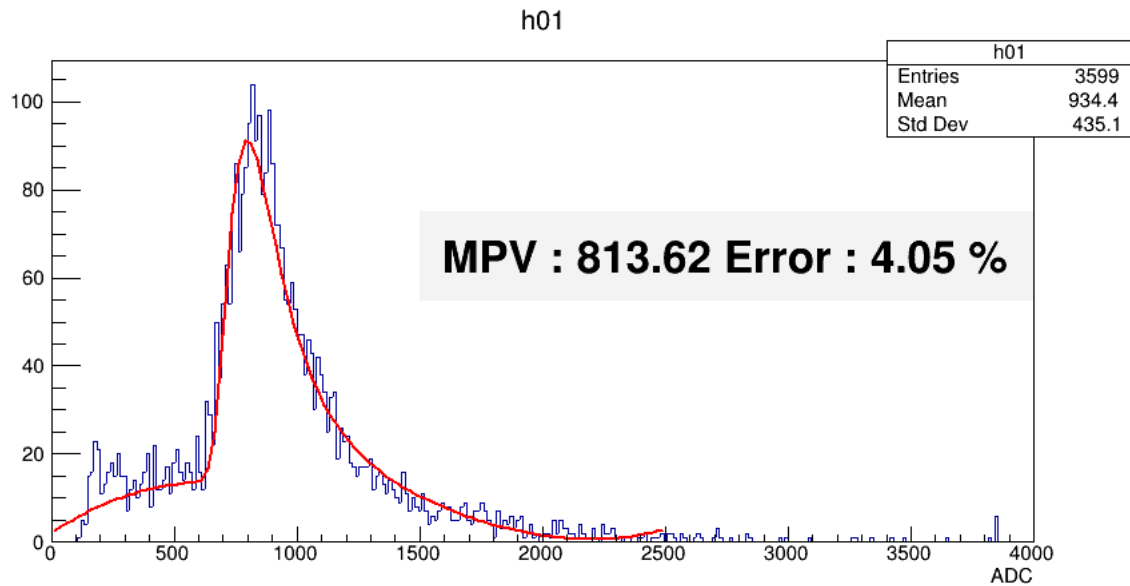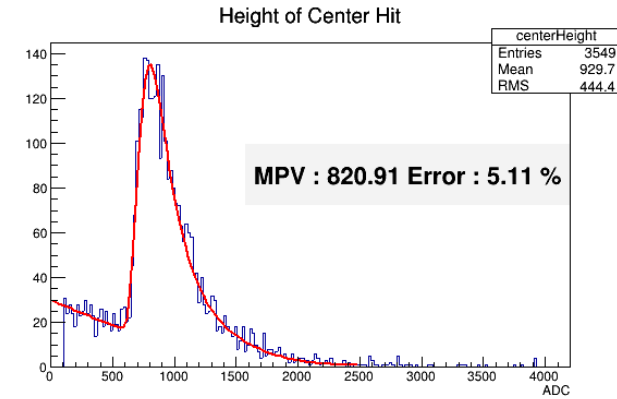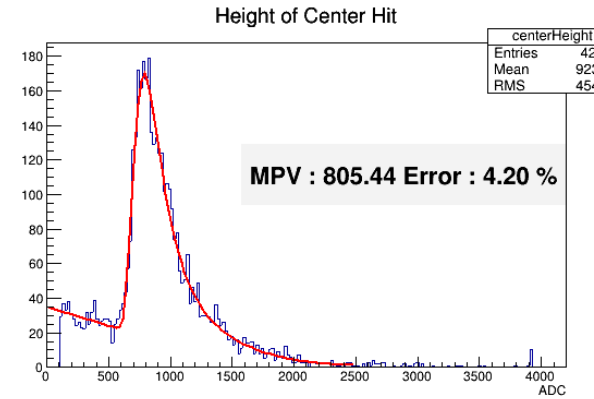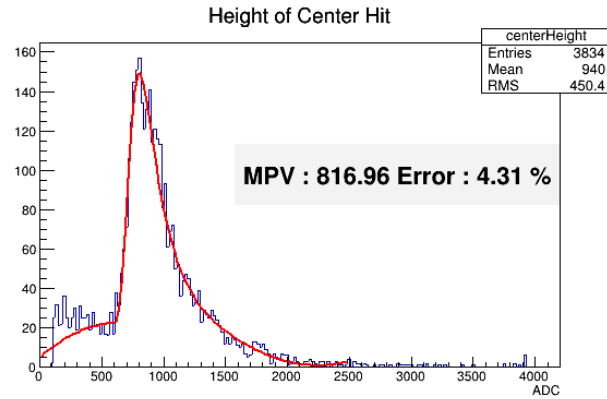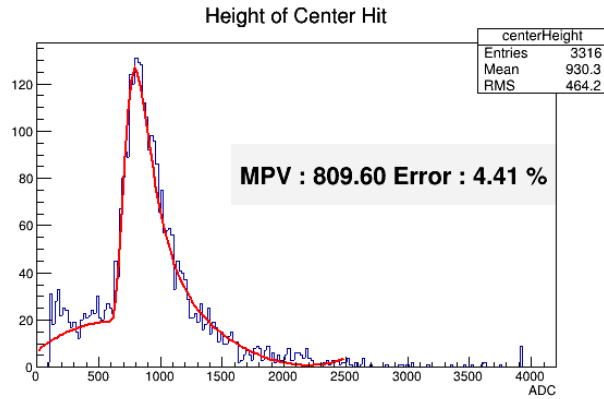



Date	MPV (Count)
18/01/10	809
18/01/16	816
18/01/23	805
18/01/30	820
19/02/22	813

Bar 19 L (optical grease)

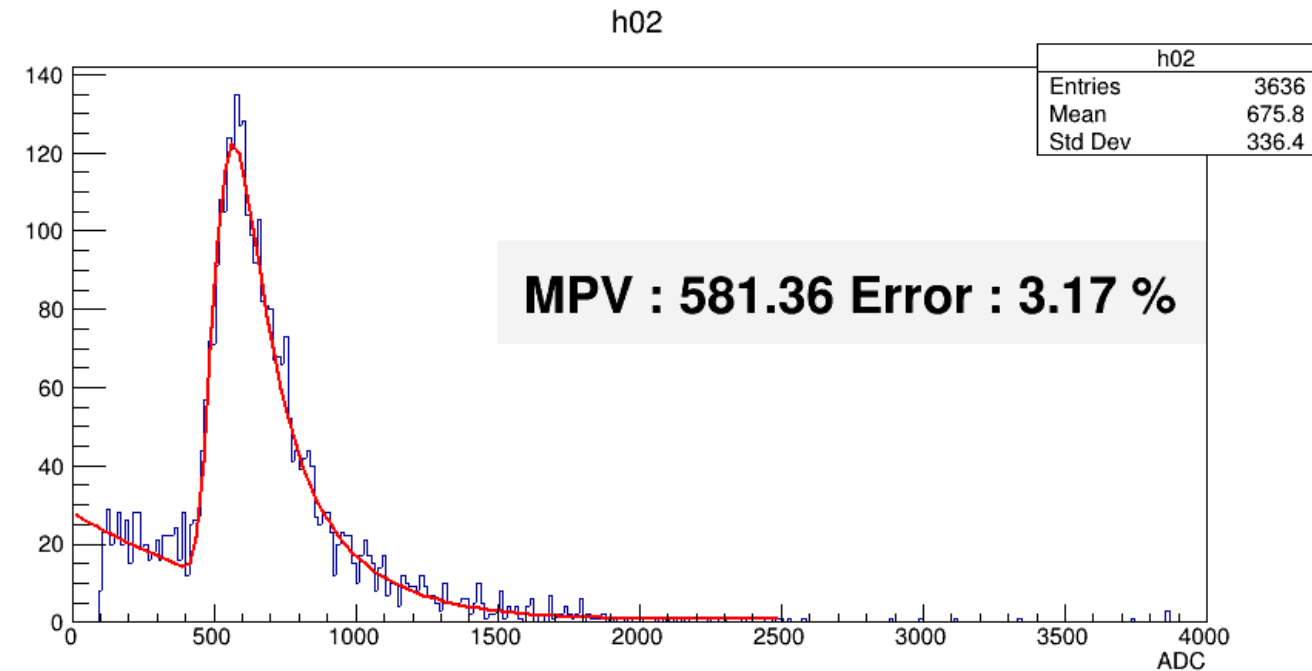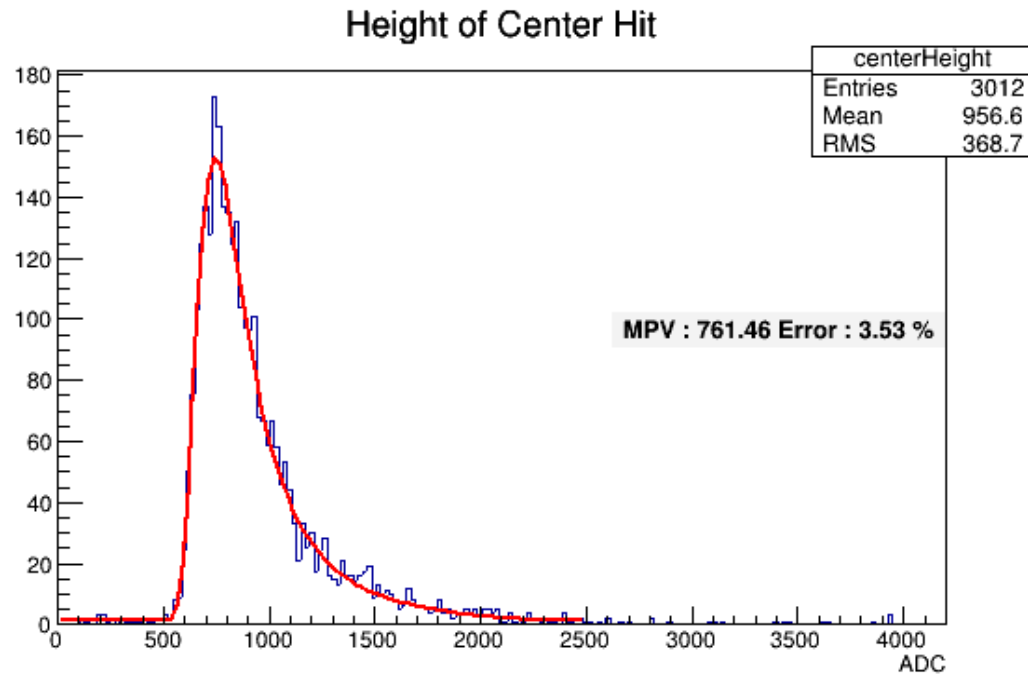

Date	MPV (Count)
18/01/3	761
19/02/22	581


Bar 19 R (optical resin)

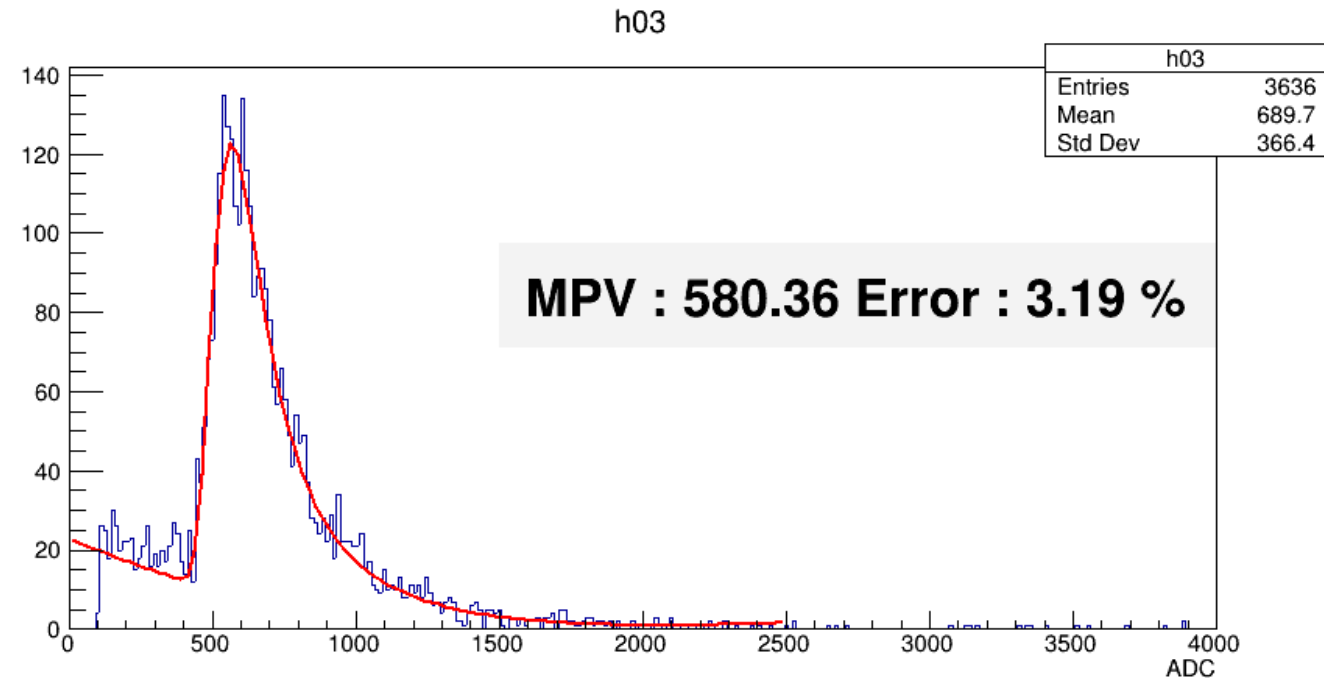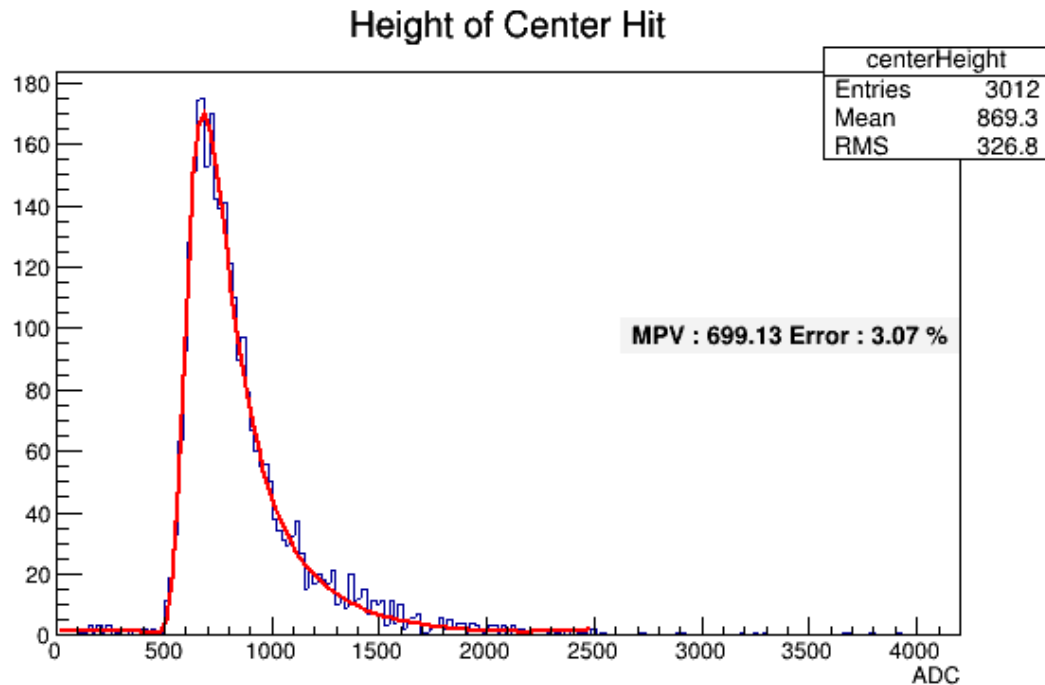

Date	MPV (Count)
18/01/3	699
19/02/22	580

Backup


Gain check(UV-3094,UV-3094-T)

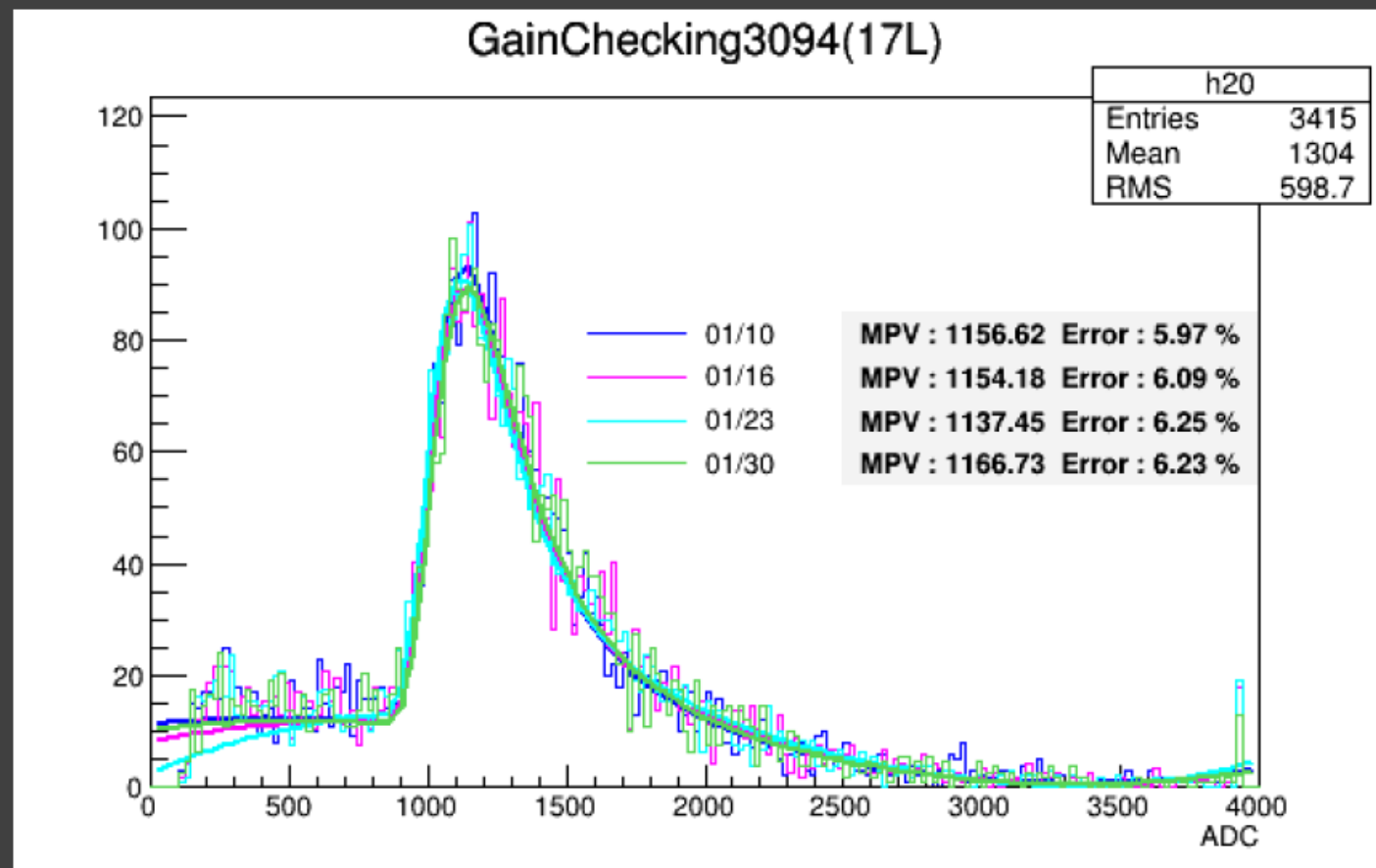

• 재진 18년 2월1일 발표 자료

Backup

Gain check(UV-3094,UV-3094-T)

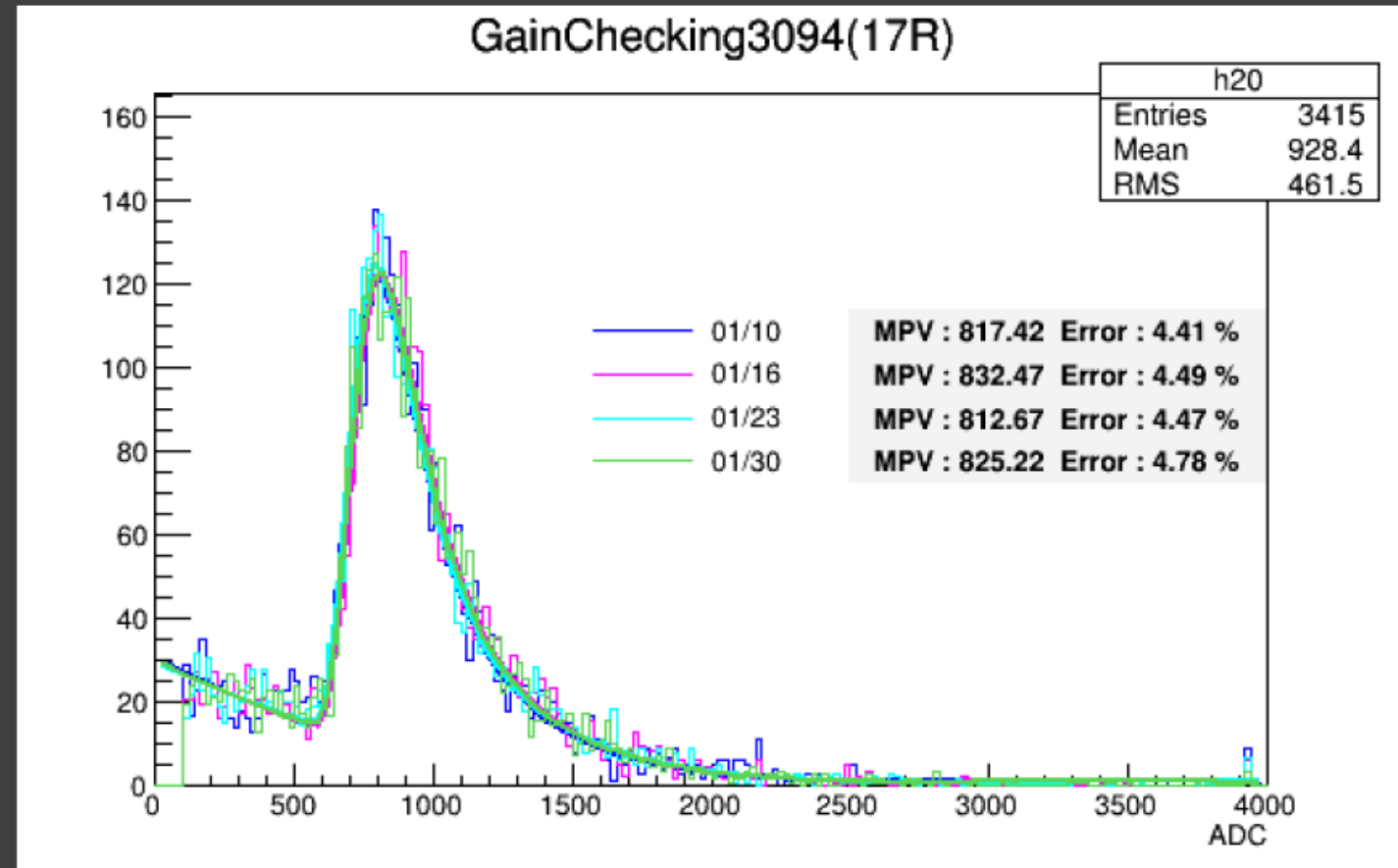

• 재진 18년 2월1일 발표 자료